



Product specification


Oilbath

ONE 10

Absolute safety for laboratory and samples, even at high temperatures through corrosion-resistant stainless steel, precise electronics and multiple temperature protection.



On this page, you can find all the essential technical data on the Memmert oil bath. Our customer relations team will be pleased to help if you want further information. If you should require a customised special solution, please contact our technical specialists at myAtmoSAFE@memmert.com.

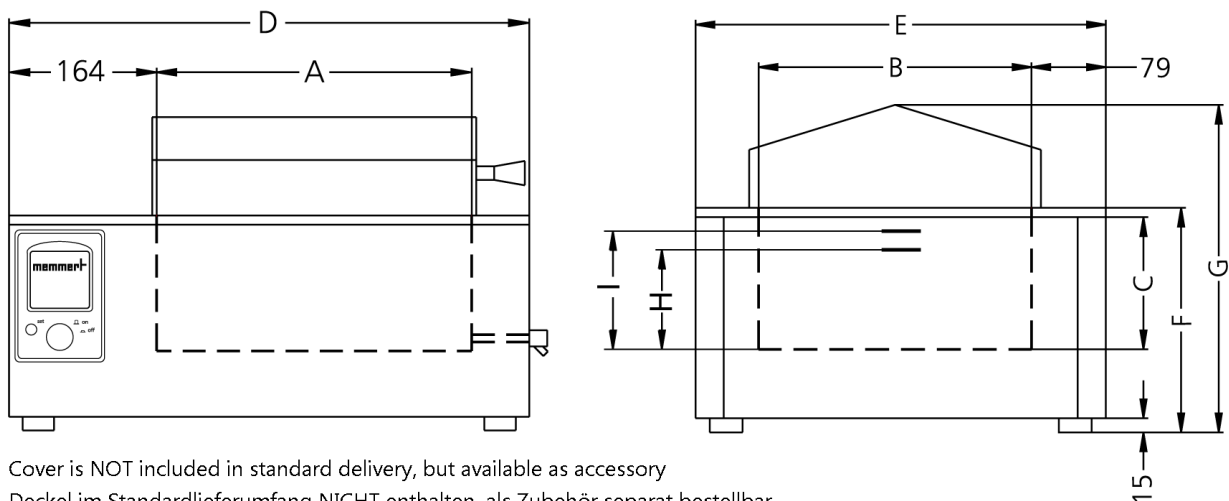

Cover is NOT included in standard delivery, but available as accessory
Deckel im Standardlieferumfang NICHT enthalten, als Zubehör separat bestellbar

Temperature

Temperature range in °C	min. 5°C above ambient up to +200°C
resolution of display and setting accuracy	0,1°C up to +99,9°C / 1°C from 100°C

Control of standard components

Controller	digital display (LED) of all set parameters, such as temperature and alarm values (0,1°C resolution) as well as time values
Timer	digital timer from 1 min. up to 999 hours for: ON (continuous operation), DELAYED ON, HOLD or HOLD set-temperature dependent with guaranteed holding time
Calibration	on controller

Safety

Temperature control	2 Pt100 sensors Class A in 4-wire-circuit, mutually monitoring and taking over the performance at the same temperature value
Temperature control	mechanical temperature limiter TB protection class 1 switching the heating off at approx. 30°C above max. temperature of the bath
Temperature control	independently working, electronic overtemperature temperature limiter TWB protection class 2
Temperature control	in case of overtemperature due to failure, the heating is switched off at approx. 10°C above the set temperature (fixed value)
Autodiagnostic system	fuzzy-supported PID microprocessor controller with integrated autodiagnostic system with fault indication
Alarm	visual and acoustic alarm at programme end and as input acknowledgement as well as in case of low liquid level, heating is switched off automatically

Heating concept

Heating Baths	corrosion-proof large-area heating on three sides
----------------------	---

Stainless steel interior

Dimensions W x H x D in mm	$w_{(A)} \times h_{(C)} \times d_{(B)}$: 350 x 140 x 210 mm
Interior	easy-to-clean interior, made of stainless steel, reinforced by deep drawn ribbing, material no. 1.4301 (ASTM 304), laser-welded
Volume	10 l
Liquid level min.	(H) 97 mm
Liquid level max.	(l) 120 mm

Textured stainless steel casing

Dimensions $w_{(D)} \times h_{(G)} \times d_{(E)}$: 578 x 337 x 356 mm

Electrical data

Voltage 230 V, 50/60 Hz

Electrical load approx. 1200 W (during heating)

Ambient conditions

Installation The vent openings in the left and back side must remain unobstructed. Minimum wall spacing on all sides is 80 mm. The minimum spacing from the top of the bath to the next ceiling is 750 mm.

Ambient temperature +5 °C to +40 °C

Humidity rh max. 80 %, non-condensing

Overvoltage category II

Pollution degree 2

Packing/shipping data

Transport information The appliances must be transported upright

Customs tariff number 8419 8998

Country of origin Federal Republic of Germany

WEEE-Reg.-No. DE 66812464

Dimensions approx incl. carton w x h x d: 670 x 450 x 380 mm

Net weight approx. 13 kg

Gross weight carton approx. 19 kg

Standard units are safety-approved and bear the test marks

